



Zimmer®  
Pulsavac® Plus AC  
Wound  
Debridement  
System



Superior care for superior debridement

# Pulsavac® Plus AC...

## Let The Power Flow Through

The Pulsavac® Plus AC Wound Debridement System combines convenience and power for a smooth and efficient flow of pulsatile irrigation and debridement. Orthopaedic surgeons have come to depend on the Pulsavac Plus System for superior debridement and wound care. Now, surgeons can utilize the benefits of the Pulsavac Plus System with the new Pulsavac Plus AC Wound Debridement System.

The Pulsavac Plus AC System is a disposable, environmentally friendly system that includes a lightweight hand piece; a variety of tips; and a reusable AC adapter box, available in nine different country configurations, that eliminates the need for alkaline batteries. Trigger-controlled pressure settings offer powerful cleaning action for bone or a gentle lavage for soft tissue. Concurrent suction and irrigation thoroughly remove debris and contaminants without flooding the field. The Pulsavac Plus AC System offers the versatility and reliability needed for effective, thorough, and environmentally friendly debridement.


### The Feel Of Total Control

#### Wide Selection of Debridement Tips

The Pulsavac Plus AC System includes a variety of debridement tips for orthopaedic applications. These tips include the dual-cone splash shield, designed to minimize splashback and help control the dangers associated with transient spray. All Pulsavac Plus AC Tips offer concurrent concentric suction to quickly remove irrigant and contaminated fluids from the operative site, providing a highly effective method of debridement and waste removal.

#### The Complete Package

The Pulsavac Plus AC System is available in pre-assembled, sterile kits for specific applications and plug configurations (see ordering information on back panel).

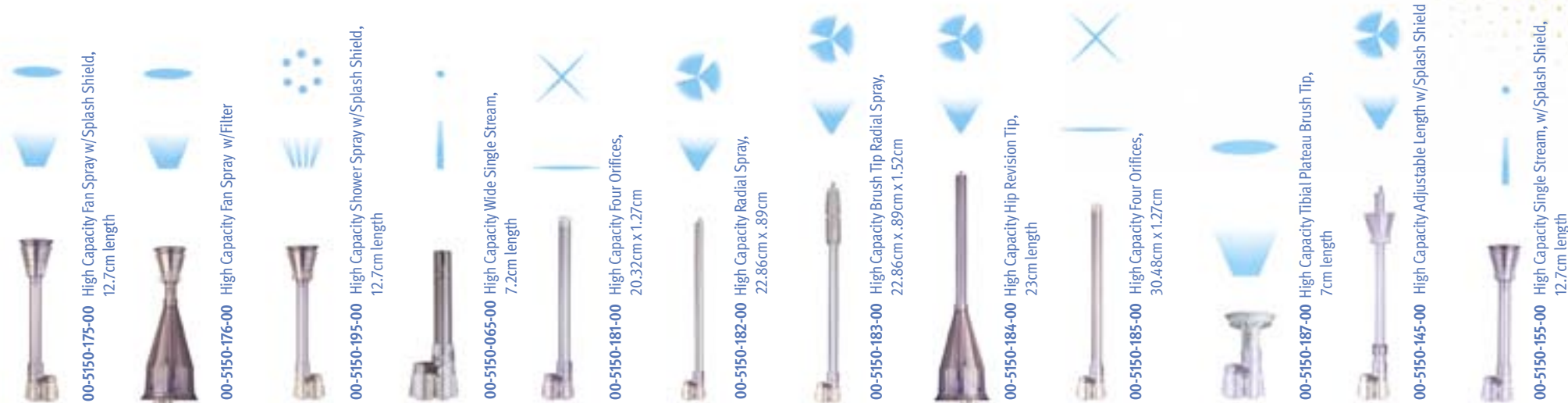

**AC Adapter**  
eliminates the need for alkaline batteries, reducing the amount of waste.



AC

**Tip Retaining Mechanism**  
provides a positive lock for secure attachment of tips.

**Textured Surface**  
near tip to help provide a sure grip and firm control, even in awkward positions.

**Rocker Switch**  
provides convenient pressure control at your fingertips.

**Textured Pistol-Grip Handle**  
provides a secure, comfortable grip to facilitate control.



### Orthopaedics

High power efficiently clears bone and cement debris during total joint procedures. Specific tips are available for hip and knee applications.

**Ergonomic, Lightweight Design**  
fits comfortably into the palm of your hand and facilitates use in tight areas.

**Independent Fluid and Suction Tubes**  
provide concentric fluid delivery and suction.

# Specifications

## Kit Packs (AC Power Pack not included\*) (Box 1)

Prod. No.	Description
00-5150-426-01	Component Kit
00-5150-476-01	Fan Spray Kit
00-5150-486-01	Hip Kit
00-5150-496-01	Shower Spray Kit

## Multi-Orifice Tips (Box 10)

Prod. No.	Description
00-5150-175-00	High Capacity Fan Spray Tip w/Splash Shield, 12.7cm length
00-5150-375-00	High Capacity Fan Spray Tip w/o Splash Shield
00-5150-176-00	High Capacity Fan Spray Tip w/Filter
00-5150-195-00	High Capacity Shower Spray Tip w/Splash Shield, 12.7cm length
00-5150-395-00	High Capacity Shower Spray Tip w/o Splash Shield

## Debridement/Suction Tips (Box 10)

Prod. No.	Description
00-5150-155-00	High Capacity Single Stream Tip w/Splash Shield, 12.7cm length
00-5150-065-00	High Capacity Large Fluid Orifice Tip w/Wide Single Stream, 7.2cm length

## Femoral Tips (Box 10)

Prod. No.	Description
00-5150-181-00	High Capacity IM Tip Four Orifices, 20.32cm x 1.27cm
00-5150-182-00	High Capacity IM Tip w/Radial Spray, 22.86cm x .89cm
00-5150-183-00	High Capacity IM Brush Tip w/Radial Spray, 22.86cm x .89cm x 1.52cm
00-5150-184-00	High Capacity IM Hip Revision Tip, 23cm length
00-5150-185-00	High Capacity IM Tip Four Orifices, 30.48cm x 1.27cm
00-5150-187-00	High Capacity Tibial Plateau Brush Tip, 7cm length


## Specialty Tip (Box 10)

Prod. No.	Description
00-5150-145-00	High Capacity Adjustable Length Tip w/Splash Shield

## Miscellaneous (Box 10)

Prod. No.	Description
00-5150-188-00	18cm Splash Shield
00-5150-189-00	23cm Splash Shield

## AC Power Pack (Reusable) (Box 1)

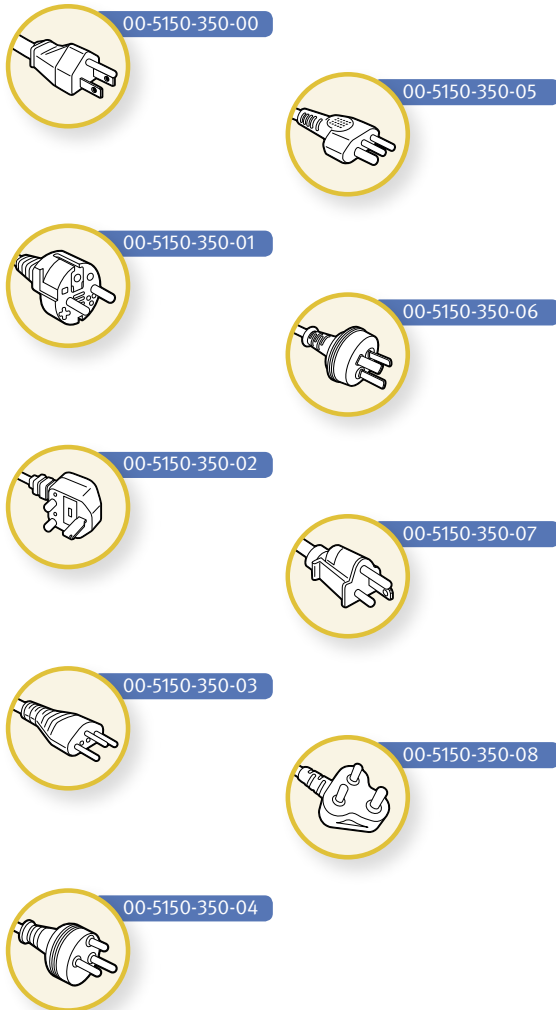

Contact your Zimmer representative or visit us at [www.zimmer.com](http://www.zimmer.com)

\*Country specific AC Power Pack sold separately. Please specify appropriate AC Power Pack with purchase of *Pulsavac* Plus AC Kit and Tips.

